

MICROCLIP/ SUPERFINE MICROCLIP

GEM
Global Excellence in Microsurgery



The right size for the Microsurgeon

Our line of hemostatic MICROCLIPS ensures a secure closure of even the most delicate vessels. These titanium clips have a unique chevron shape and interlocking atraumatic teeth that prevent multi-directional slippage.

Pure Titanium Annealed to a Malleable Consistency

- Provides maximum vessel occlusion
- Non-ferrous for use with MRI or related scanning equipment
- Convenient adhesive backing on each cartridge




Rx ONLY. For safe and proper use of the device, please refer to full Instructions for Use.

MICROCLIP/ SUPERFINE MICROCLIP



Chevron Shape

Provides progressive and encompassing tip-to-tip closure before complete and secure occlusion of tissue or vessels

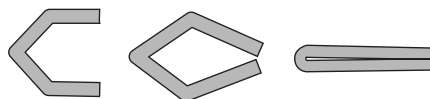


Diamond Shape Grooving

- Diamond grooving and interlocking atraumatic teeth design of the clip ensures enhanced security in vessel closure
- Striation is composed of small pyramidal heads, similar to tungsten carbide inserts in needle holders, ensuring a firm multi-directional non-crushing grip



Measurements

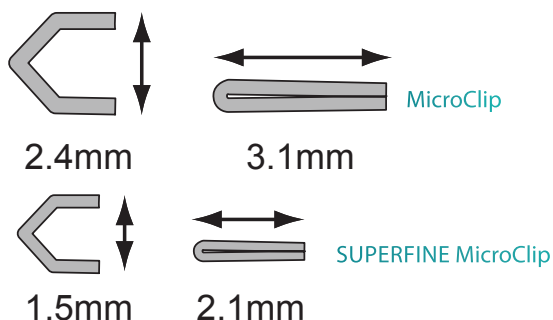

Ordering Information

MicroClips	GEM2431	Box of 30
MicroClip Appliers – 15 cm	GEM615-1G	Each
MicroClip Appliers – 20 cm	GEM620-1G	Each
SUPERFINE MicroClips	GEM1521	Box of 30
SUPERFINE Appliers – 15 cm	GEM715-SF	Each
SUPERFINE Appliers – 20 cm	GEM720-SF	Each

Manufactured by Péters Surgical and Distributed in the USA by Synovis Micro Companies Alliance, Inc.

Synovis Micro Companies Alliance

GEM design and Superfine are trademarks of Baxter International Inc., or its subsidiaries.

FAX: 205-341-3192 | EMAIL: scma_information@baxter.com | PHONE: 800-510-3318 | 439 Industrial Lane, Birmingham, AL 35211
Synovis Micro Companies Alliance, Inc. (a subsidiary of Baxter International Inc.)

USMP/MG219/18-0007